

NEW

Electronically rotating thumb glides effortlessly into grip position

18 different grip options available




touch
bionics BY ÖSSUR

i-limb® ultra

An externally powered, multi-articulating prosthetic hand with five individually powered digits and powered thumb rotation, with manual override.

- Individually motorized digits with compliant grip and stall-out ability
- Speed boost:
Adjustable open/close speed for digits
- Auto-grasp:
Anti drop safety feature for added reassurance
- Vari-grip:
Feature for enhanced grip strength
- Multiple covering options available to suit individual lifestyles
- Titanium digits
Check out our titanium strengthened digits which increase the maximum carry load by 50%. Available in S/M/L.

METHODS OF CONTROL:

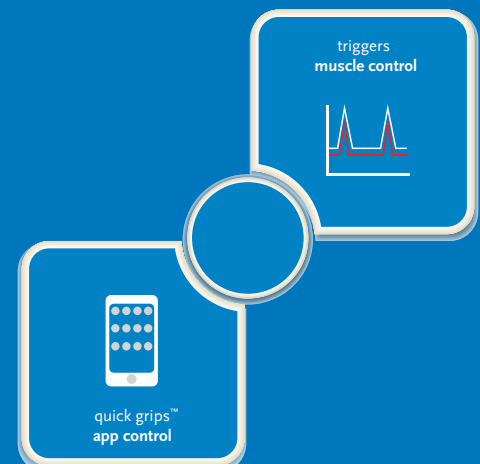



ÖSSUR

TOUCH
SOLUTIONS

Össur Touch Solutions are designed to ensure the best possible outcomes for people with upper-limb deficiencies.

Technical Information

Voltage	7.4 V (nominal)			
Max. Current	5 A			
Battery Capacity	Rechargeable lithium polymer; 7.4 V (nominal); 2000 mAh capacity; 1300 mAh capacity			
Max. hand load limit (static limit)	40kg/88lb (Extra Small) 90kg/198lb (Small/Medium/Large)			
Finger carry load (static limit)	20kg/44lbs (Extra Small) 32kg/71lbs (Small/Medium/Large)			
Time from open position to full power grip	0.8 seconds			
Device Weight		Extra Small	Small	Medium/Large
	QWD	472g/1.04lb	512g/1.13lbs	528g/1.16lbs
	WD	432g/0.96lbs	472g/1.04lbs	488g/1.08lbs
	Flexion Wrist	572g/1.26lbs	612g/1.35lbs	628g/1.38lbs
	Friction Wrist	467g/1.03lbs	507g/1.12lbs	523g/1.15lbs

Note: Titanium digits add an additional 30g/0.07lbs per hand

Part #	Description
PL403035	QWD, Large, RH
PL403034	QWD, Large, LH
PL403033	QWD, Medium, RH
PL403032	QWD, Medium, LH
PL403031	QWD, Small, RH
PL403030	QWD, Small, LH
PL403029	QWD, Extra Small, RH
PL403028	QWD, Extra Small, LH

Part #	Description
PL403007	Wrist Disarticulation, Large, RH
PL403006	Wrist Disarticulation, Large, LH
PL403005	Wrist Disarticulation, Med., RH
PL403004	Wrist Disarticulation, Med., LH
PL403003	Wrist Disarticulation, Small, RH
PL403002	Wrist Disarticulation, Small, LH
PL403001	Wrist Disarticulation, XSmall, RH
PL403000	Wrist Disarticulation, XSmall, LH

Part #	Description
PL403047	Flexion Wrist, Large, RH
PL403046	Flexion Wrist, Large, LH
PL403045	Flexion Wrist, Medium, RH
PL403044	Flexion Wrist, Medium, LH
PL403043	Flexion Wrist, Small, RH
PL403042	Flexion Wrist, Small, LH
PL403049	Flexion Wrist, Extra Small, RH
PL403048	Flexion Wrist, Extra Small, LH

FOLLOW ÖSSUR ON



USA (800) 233-6263
WWW.OSSUR.COM

CANADA (800) 663-5982
WWW.OSSUR.CA

