

ICEROSS SEAL-IN® V TRANSTIBIAL LINER WITH WAVE




Iceross Seal-In V features the latest Seal-In membrane for the best in-socket seal technology available today.

Based on the Iceross Cushion profile, Iceross Seal-In V has a soft distal end for close fit and comfort.

Iceross Seal-In V is available with the seal in two different locations, Standard and High Profile. The High Profile variant should be used when residual limb length allows for it, especially in case of sensitive spots on the distal tibial crest.

Össur recommends Iceross Seal-In V for amputees already past the initial phase of volume and shape conditioning, during the first months post-amputation Iceross Dermo Cushion & Iceross Sleeve are recommended.

ICEROSS®
Trans tibial – Seal-In®

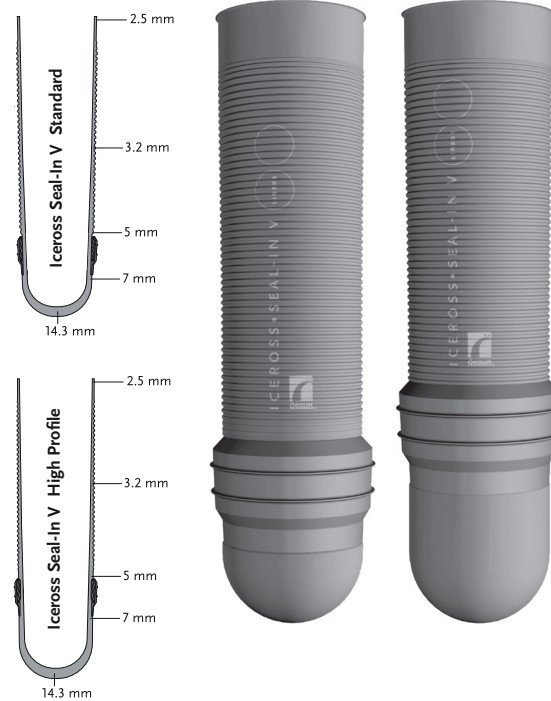

USER INFORMATION	
Amputation Level:	Trans tibial
Impact Level:	All
LINER INFORMATION	
Sizes:	18,20,22,23.5,25,26.5,28,30,32,34,36.
Profile:	3mm wave

SUSPENSION METHOD



SEAL-IN V


FEATURES



FABRIC COVER ACTIVE SKIN CARE WAVE SILKEN INNER SURFACE

SEAL HEIGHTS AND TRIM LINES

LINER SIZE	MINIMAL RESIDUAL LENGTH*	
	STANDARD	HIGH PROFILE
18-20	10cm	13cm
21-25	11cm	14cm
26-30	12cm	15cm
31-36	13cm	16cm

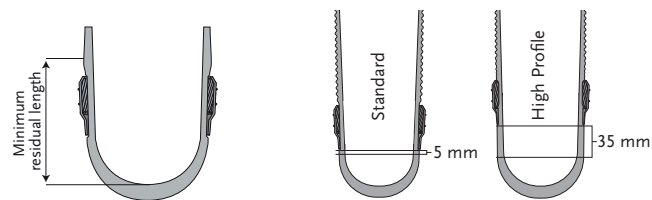

*Residual limb shapes and trim line preference varies, minimum residual length values are approximates only.

PART NUMBER INFORMATION - ICEROSS SEAL-IN V TRANSTIBIAL LINER WITH WAVE

PART#	PROFILE	DESCRIPTION
I-4713XX	Wave 3mm	Seal-In V
I-4723XX	Wave 3mm	Seal-In V High Profile

www.ossur.com



Life Without Limitations®