



FLEX-FOOT® ASSURE + UNITY™ SLEEVELESS VACUUM

More than just a comfortable foot, the easy-to-fit Flex-Foot Assure is a simple yet dynamic solution for patients who need the benefits of performance in an easy-to-walk design. The ultra lightweight, cushioned heel reduces impact, while the full-width torsional keel optimizes walking symmetry and protects the integrity of the sound limb.



And now Flex-Foot Assure users can go sleeveless, with our Unity Sleeveless Vacuum System, for moisture evacuation and volume control.

Ask your Össur rep for a demo today, or visit ossur.com/low-activity to schedule an on-site course and earn 4.0 CEUs.



FLEX-FOOT® PRODUCT SELECTOR

<< LOWER IMPACT | HIGHER IMPACT >>



IMPACT LEVELS



FEATURES



ACTIVE TIBIAL PROGRESSION

Vertical forces generated at heel contact are initially absorbed and translated into forward motion described as Active Tibial Progression.



FULL LENGTH TOE LEVER

The full-length keel/toe lever matches the length of the sound foot, producing a smoother and more symmetrical gait, while maximizing energy return.




PROPORTIONAL RESPONSE

The layering of carbon fiber optimizes the deflection of the forefoot from mid-stance to toe-off proportionally to the user's weight and impact level.

UNITY FOR FLEX-FOOT® ASSURE

SUGGESTED BILLING CODE*
L-5981, L-5781

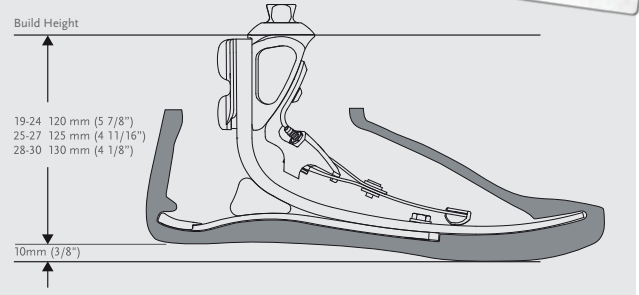

RATINGS

- Low to moderate impact levels (community ambulator)
- Patient Weight: 99-300 lbs (45-136 kg)

OPTIONS

- Available in sizes 19-30, categories 1-5
- Male pyramid or 30mm tube
- FAPUXYYZ* Flex-Foot Assure with Unity
- FAXEYXXZ* Flex-Foot Assure
- FAPEYXXZ* Flex-Foot Assure with pyramid
- Foot comes with foot cover

UNITY KITS

- UFA0024 FAP-JBP Unity Pump S22-24
- UFA0027 FAP-JBP Unity Pump S25-27
- UFA0030 FAP-JBP Unity Pump S28-30
- UTT0001 Unity TT Valve Kit
- UTF0001 Unity TF Valve Kit

24-MONTH WARRANTY ON UNITY PUMP KIT

12-MONTH WARRANTY ON UNITY VALVE & MEMBRANE

36-MONTH WARRANTY ON FLEX-FOOT ASSURE

*X: category (1-5); YY: size (19-30); Z: side (L/R)
Add "BR" to the part # for a brown foot cover

CATEGORY SELECTION GUIDE FOR FLEX-FOOT ASSURE

WEIGHT	FOOT SIZE											
	19	20	21	22	23	24	25	26	27	28	29	30
45 - 61kg (99 - 135lbs)	1	1	1	1	1	1	1	1	1	1		
62 - 80kg (136 - 176lbs)	2	2	2	2	2	2	2	2	2	2	2	2
81 - 97kg (177 - 214lbs)					3	3	3	3	3	3	3	3
98 - 114kg (215 - 251lbs)					4	4	4	4	4	4	4	4
115 - 136kg (252 - 300lbs)							5	5	5	5	5	5

* The responsibility for accurate coding lies solely with the provider treating the patient. Össur assumes no responsibility or liability for the provider's coding decisions. Össur's coding suggestions rest on its best judgment and are subject to revision.

FOLLOW ÖSSUR ON



TEL (800) 233-6263
FAX (800) 831-3160
WWW.OSSUR.COM

