



ÖSSUR[®]

LIFE WITHOUT LIMITATIONS

REBOUND[®] DUAL

Versatility and
stability delivered





Introducing Rebound® DUAL

Össur's Rebound DUAL knee brace delivers functional support for ligament instabilities, including ligament instability with osteoarthritis, in a single, versatile brace. Suitable for low-to-high impact activities as well as sliding sports, Rebound DUAL is your new, go-to knee brace.

Slide-to-size

The multi-purpose upper frame can be adjusted from 12.5" to 14" to suit the patient's height and help ensure the perfect fit. It can also be adjusted for varus/valgus alignment.



Comfort plus

The lightweight aluminum frame allows for easy adjustability by the clinician, as well as optimal comfort for the patient. Doeskin thigh liners and easy-on, easy-off strapping are designed to increase patient satisfaction.

No-slip ActiveGrip™

The breathable ActiveGrip calf liner has anti-migration capabilities. Once in place, it stays in place. Even the most active users can go all out with confidence.



Adjustable SmartDosing®

The optional SmartDosing system offers dynamic tension controls, allowing fine-tuned adjustments for osteoarthritis pain relief as needed by the user.



REBOUND® DUAL

Flexible, comfortable and reliable support for improved ligament stability.

Adjustable frames means fewer parts to order

- Pre-fabricated version comes in 7 sizes and right or left configurations
- Custom options also available

FEATURES AND BENEFITS

- 1 Slide-to-size adjustable upper frame
Adjustable from 12.5"-14"
- 2 Cold-malleable aluminum frame
Easy to contour to patient's unique anatomy
- 3 Simple strapping system
Easy-on, easy-off for the patient
- 4 Field-serviceable D-rings
Can be adjusted to reduce pressure on the patient
- 5 ActiveGrip™ calf liner
Breathable with anti-migration capabilities
- 6 Doeskin thigh liner
Increased comfort on the thigh
- 7 Cut-to-fit strap pads
Trim to optimize patient comfort

REBOUND® DUAL ORDERING INFORMATION

PART#	SIZE	MEASURING ^{1,2}
B-242X10001	XS	Caliper: 3"-3.5" (7.6-8.9 cm) Circ: 15"-16.5" (38.1-41.9 cm)
B-242X10002	SMALL	Caliper: 3.5"-4" (8.9-10.2 cm) Circ: 16.5"-18.5" (41.9-47 cm)
B-242X10003	MEDIUM	Caliper: 4"-4.5" (10.2-11.4 cm) Circ: 18.5"-20.5" (47-52.1 cm)
B-242X10004	LARGE	Caliper: 4.5"-5" (11.4-12.7 cm) Circ: 20.5"-22.5" (52.1-57.2 cm)
B-242X10005	XL	Caliper: 5"-5.5" (12.7-14 cm) Circ: 22.5"-24.5" (57.2-62.2 cm)
B-242X10006	XXL	Caliper: 5.5"-6" (14-15.2 cm) Circ: 24.5"-27.5" (62.2-69.9 cm)
B-242X10007	3XL	Caliper: 6"-6.5" (15.2-16.5 cm) Circ: 27.5"-30" (69.9-76 cm)
CUSTOM PART#	MEASURING ³	
B-155X10001	SM, DMS, CCS, CAD, CAST	
ACCESSORY PART#	MEASURING	
B-701090060	Right medial/left lateral SmartDosing Kit	
B-701090061	Left medial/right lateral SmartDosing Kit	

X = 5 (right) or 6 (left)

1 Caliper: M/L at joint width in a standard, weight-bearing position

2 Thigh Circumference: 6" above mid-patella






3 SM = SmartMeasure, DMS = Digital Measuring System, CCS = Coordinate Cast System, CAD = Computer Aided Design




CUSTOM PAINT COLORS:

-  Black (default)
-  Electric Blue
-  Metallic Silver
-  White
-  Red

CUSTOM LINER COLORS:

-  Metallic Silver (default)
-  Electric Blue
-  Yellow
-  Purple
-  Black

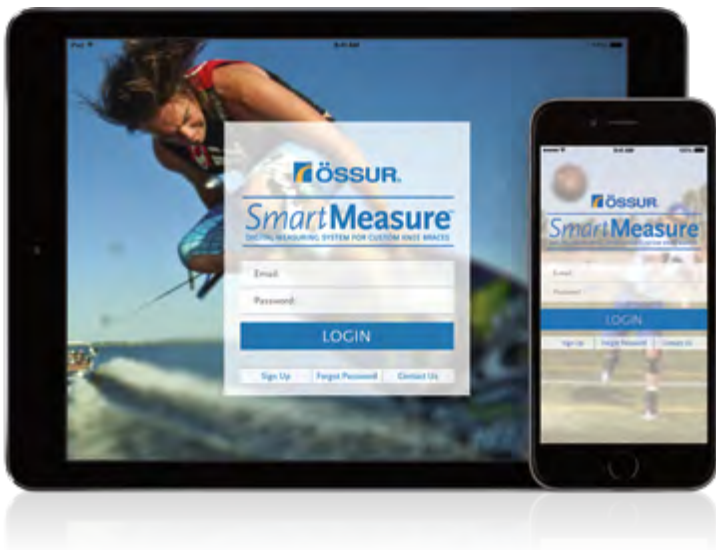

Nothing measures like SmartMeasure™

An app for iPhone® & iPad® designed to order pre-fabricated knee braces or measure custom knee braces. SmartMeasure is Össur's most intuitive, efficient and accurate way to measure and order Össur knee braces, including Rebound DUAL, Rebound PCL and Rebound Cartilage.



Apple, the Apple logo, iPhone and iPad are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc.

FOLLOW ÖSSUR ON



TEL (800) 233-6263
CANADA (800) 663-5982
WWW.OSSUR.COM

