

REBOUND® DUAL

The Rebound DUAL knee brace delivers functional support for ligament instabilities, including ligament instability with osteoarthritis, in a single, versatile aluminum brace. Suitable for low-to-high impact activities as well as sliding sports, Rebound DUAL is your new, go-to knee brace.



INDICATIONS FOR USE

- Supports ACL, MCL, LCL, PCL or combined instabilities
- Ligament instabilities with osteoarthritis


FEATURES AND BENEFITS

- Polycentric hinge for stability, durability, and low-profile design
- Adjustable tibia strap captures tibia for optimal stability
- Optional modular SmartDosing® kit available when needed for fine-tuned pain management and offloading
- Breathable ActiveGrip™ calf liner with anti-migration capabilities

REQUIRED MEASUREMENTS

- Caliper measurement: Medio-laterally at joint width in a standard, weight-bearing position
- Circumference: Just below the distal border of the patella

CUSTOM BRACE COLORS

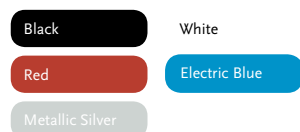

CUSTOM LINER



REBOUND® DUAL

ORDERING INFORMATION

REBOUND DUAL

Size	Caliper	Circumference	Right	Left
XSmall	3-3.5" (7.6-8.9 cm)	15-16.5" (38.1-41.9 cm)	B-242510001	B-242610001
Small	3.5-4" (8.9-10.2 cm)	16.5-18.5" (41.9-47 cm)	B-242510002	B-242610002
Medium	4-4.5" (10.2-11.4 cm)	18.5-20.5" (47-52.1 cm)	B-242510003	B-242610003
Large	5-5.5" (12.7-14 cm)	22.5-24.5" (57.2-62.2 cm)	B-242510004	B-242610004
X-Large	5.5-6" (14-15.2 cm)	24.5-27.5" (62.2-69.9 cm)	B-242510005	B-242610005
XX-Large	6-6.5" (15.2-16.5 cm)	27.5-30" (69.9-76 cm)	B-242510006	B-242610006
XXX-Large	6-6.5" (15.2 – 16.5 cm)	27.5-30" (69.9 – 76 cm)	B-242510007	B-242610007

M/L at joint width in a standard, weight-bearing position 2 Thigh Circumference: 6" above mid-patella

REBOUND DUAL CUSTOM

Measuring	Right	Left
SmartMeasure, DMS, Coordinate Cast System, Computer Aided Design	B-155510001	B-155610001

REBOUND DUAL ACCESSORIES

Description	Part #
PCL Strap Kit XS-XL	B-701090057
PCL Strap Kit X-XXXL	B-701090058
Field Serviceable D-Rig with Rivets Kit	B-701090059
Smart Dosing Kit RM/LL	B-701090060
Smart Dosing Kit RL/LM	B-701090061



SMARTMEASURE®

A digital measuring system for custom knee braces on the iPhone or iPad. The most intuitive, efficient and accurate way to measure for any Össur custom knee brace, including the CTi. Utilizes a section-based, step-by-step approach to guide you through the entire brace selection, measuring and order process.



DIGITAL MEASUREMENT SYSTEM (DMS)

This system simply requires minimal physical measurements and two digital photos, which can then be sent to customer service. This system is a quick, easy and reliable method of measuring for a CTi without the need for training.